


## Vibratome Blades

£22.50

### GS-053234

Available in packs of 50, these single edged blades are designed to easily fit all versions of the Vibratome range, allowing sectioning down to 10-20 micron range with good technique and fixed material.

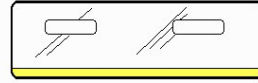

## Sapphire Knife

£740.00

### GS-053237

Developed specifically for the Vibratome, this sapphire microtomy knife makes it possible for the first time to cut thin Vibratome sections with microscopic resolution similar to that obtained with thick (1 micron) plastic sections. The sapphire knife allows serial sectioning with a Vibratome 1000 of a thickness in the 6-10 micron range with good technique and fixed material, revealing subcellular features never before seen in Vibratome sections. The sapphire knife can be resharpened.



## Foot Pedal

£135.00

### GS-053300

The foot pedal is an optional component available for Vibratome tissue sectioning systems. The pedal plugs into the back of all Vibratome instruments and is activated by a selector switch. Once activated, the foot switch is utilized for initiating a sectioning sequence and reversing the blade travel.



## Specimen Mounts 1/2" head

£ 32.50

### GS-053753

Specimen mounts are designed to allow for the rotation of the specimen in the Vibratome vice assembly for sectioning at different horizontal orientations.



## Specimen Mounting Blocks

£ 27.50

### GS-053747

The Vibratome comes standard with a square vice mounting device. These specimen blocks are utilized for mounting specimens and placing in the Vibratome holding vice. The blocks are anodized black for easy specimen viewing and can be reused. They are shipped three to a package.




## Rotating Stage Assembly

£230.00

### GS-053760

The rotating stage assembly is designed to give the researcher maximum flexibility in the mounting of specimens. It can be turned to any horizontal and vertical angle via a ball joint mounting arrangement.

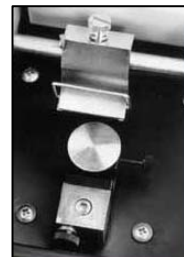

# Vibratome Accessories.

## Deluxe Specimen Bath GS-053745

£245.00

The Deluxe specimen bath mounts inside the existing black Vibratome bath. This large removable bath allows for the separation of buffer and ice kept in the original black bath. The bath can be sterilized and has round mounting disks that can be easily removed for mounting specimens. The disks and bath surface are marked with degrees of rotation for the repeatable placement and measurement of specimens.



## Section Mounting Device GS-053225

£147.50

A simple device designed to mount floated histological sections onto standard glass slides. This product facilitates the mounting of all types of thin (2-100 micron) tissue slices. It is especially useful for mounting tissue slices not supported by embedding material, such as those produced with a Vibratome, freezing microtome (cryostat) or sliding microtome. It also simplifies the mounting of tissue sections embedded in a water-miscible supporting medium such as polyethylene glycol (PEG).



## Accessories shown

|          |                              |         |
|----------|------------------------------|---------|
| GS-53747 | Specimen mounting block Pk 3 | £ 27.50 |
| GS-53760 | Rotating Stage Assembly      | £230.00 |
| GS-53766 | Rotating Mount Storage Block | £120.00 |
| GS-53750 | V-Block Adapter              | £ 22.50 |
| GS-53744 | Standard Specimen Tray       | £135.00 |



Oct 2007

R & L Slaughter Ltd. can accept orders by post, telephone, fax or e-mail. All prices are exclusive of V.A.T. and are correct at time of going to print but are subject to change without prior notice. Carriage will be charged at £10.00 (+V.A.T) per order. A low value order charge for orders less than £60.00 (excluding V.A.T) will be charged to non account holders. R & L Slaughter Ltd. accept all major credit cards.

R & L Slaughter Ltd. has been accredited to ISO 9001:2000.



R & L Slaughter Ltd. Units 11 & 12 Upminster Trading Park, Warley Street, Upminster. RM14 3PJ  
Tel 01708 227140 Fax 01708 228728 website: [www.slaughter.co.uk](http://www.slaughter.co.uk) e-mail: [sales@slaughter.co.uk](mailto:sales@slaughter.co.uk)